

PCR Kit Quality Control Certificate

Pathogen
<i>Tenacibaculum dicentrarchi</i>

Table 1: Concentration in ng/μL of the sample control and its dilutions.

	1	2	3	4	5	6	7	8	9
	1×10^{-1}	1×10^{-2}	1×10^{-3}	1×10^{-4}	1×10^{-5}	1×10^{-6}	1×10^{-7}	1×10^{-8}	1×10^{-9}
N° of copies	1×10^8	1×10^7	1×10^6	1×10^5	1×10^4	1×10^3	1×10^2	1×10^1	1×10^0

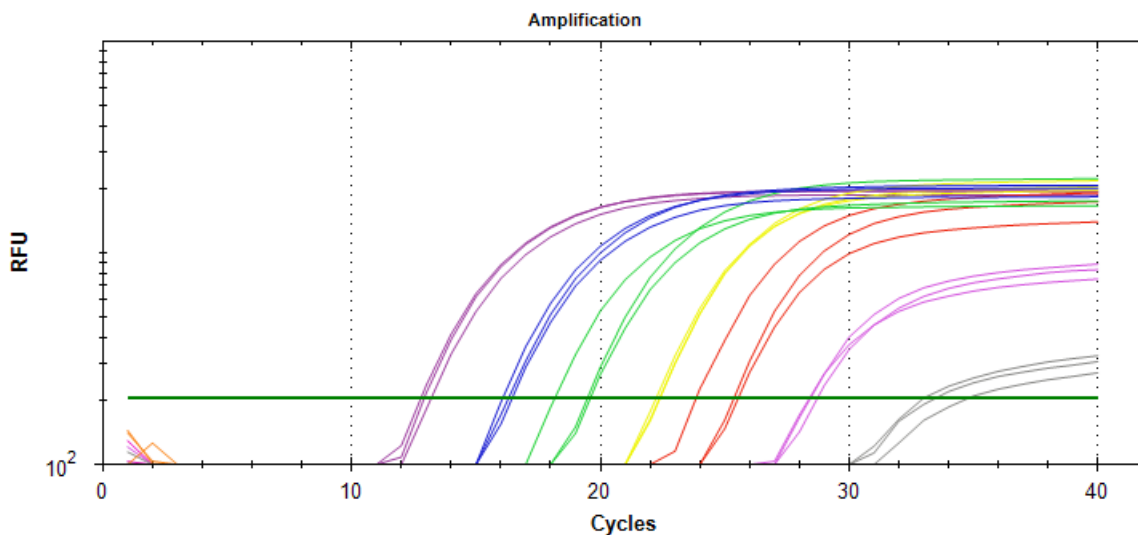


Figure 1: Logarithmic amplification curves for *Tenacibaculum dicentrarchi*

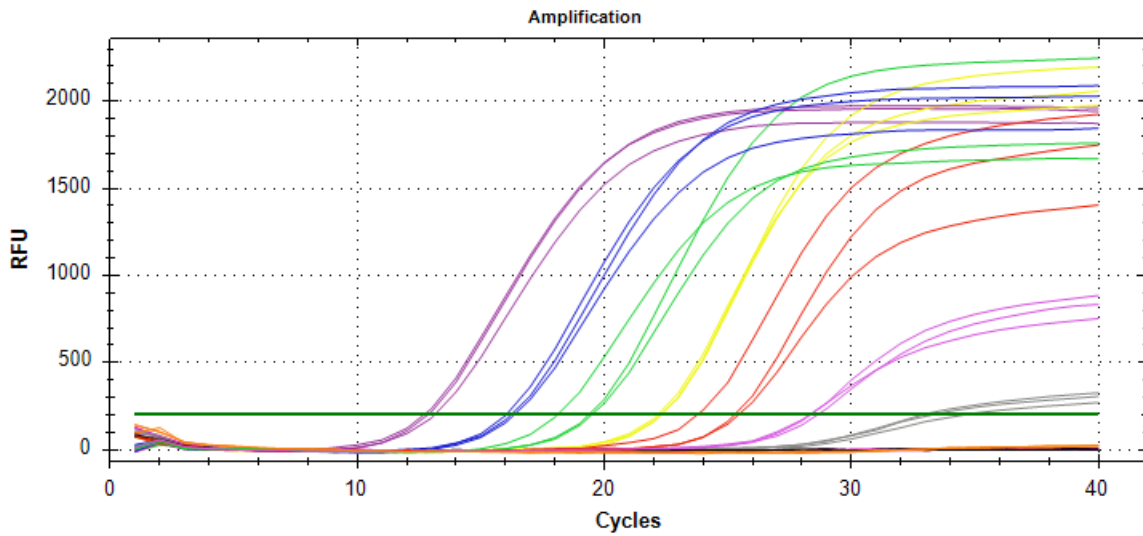


Figure 2: Standard amplification curves for *Tenacibaculum dicentrarchi*

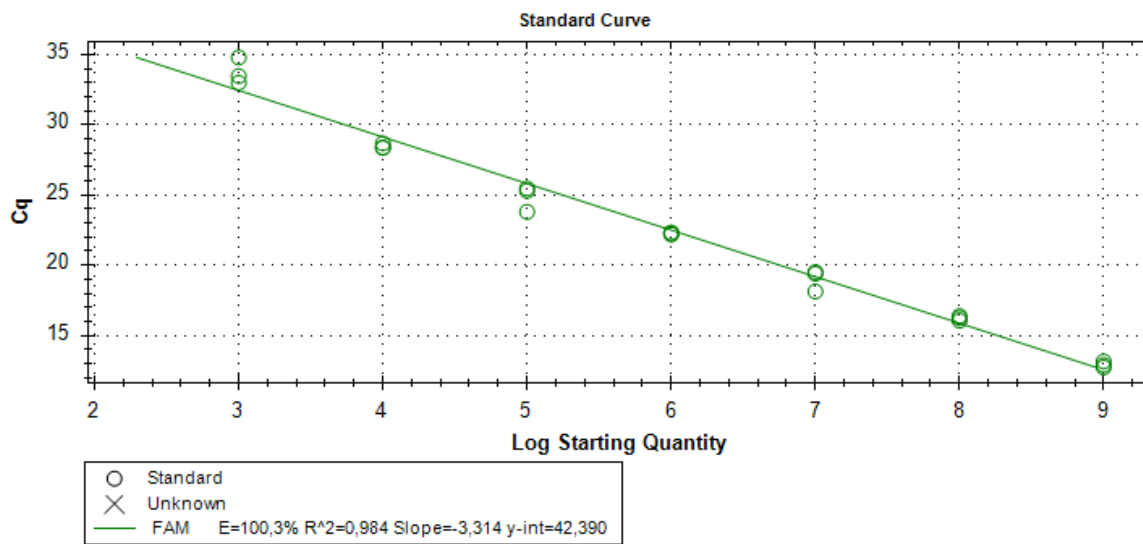


Figure 3: Standard curve for *Tenacibaculum dicentrarchi*

Pathogen
<i>Tenacibaculum maritimum</i>

Table 1: Concentration in ng/μL of the sample control and its dilutions.

	1	2	3	4	5	6	7	8	9
	1×10^{-1}	1×10^{-2}	1×10^{-3}	1×10^{-4}	1×10^{-5}	1×10^{-6}	1×10^{-7}	1×10^{-8}	1×10^{-9}
N° of copies	1×10^8	1×10^7	1×10^6	1×10^5	1×10^4	1×10^3	1×10^2	1×10^1	1×10^0

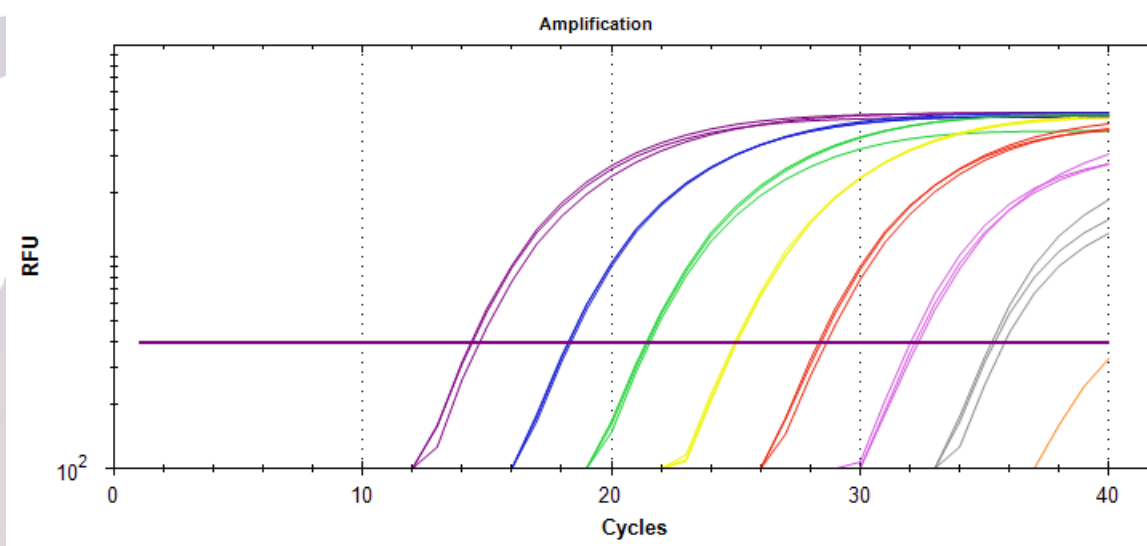


Figure 1: Logarithmic amplification curves for *Tenacibaculum maritimum*

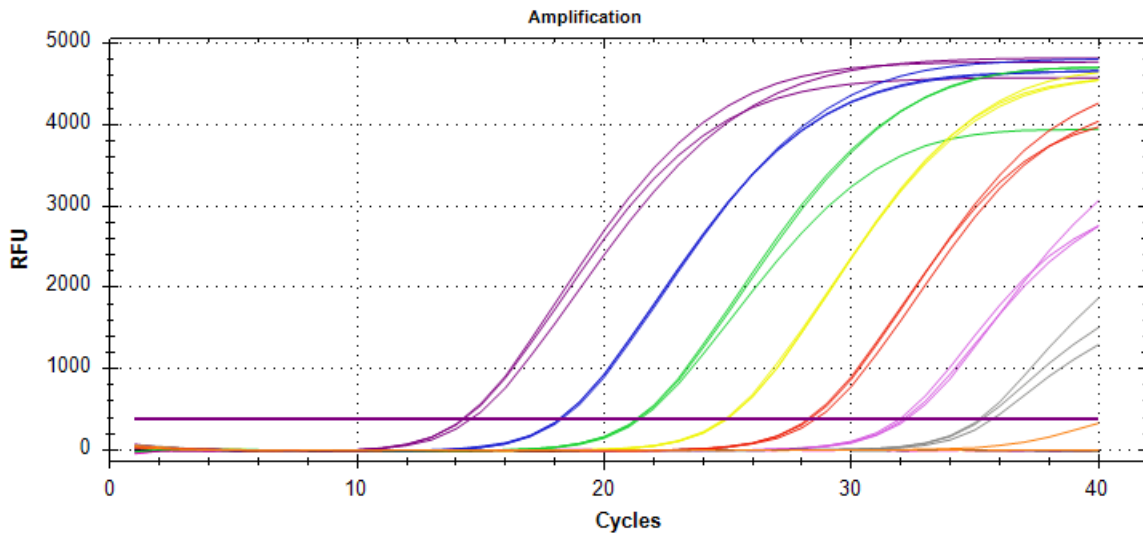


Figure 2: Standard amplification curves for *Tenacibaculum maritimum*

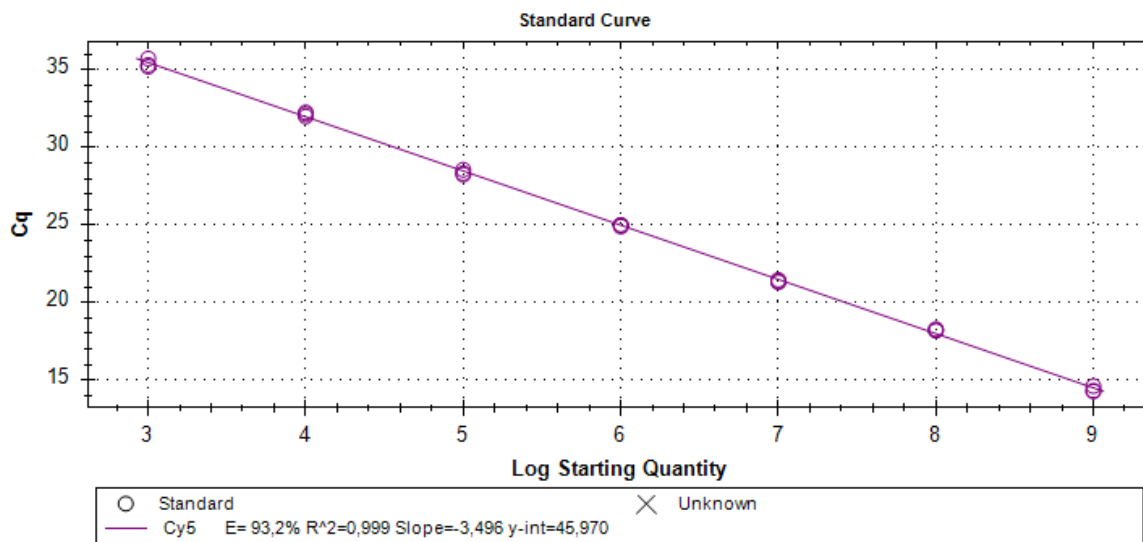


Figure 3: Standard curve for *Tenacibaculum maritimum*

Summary of Ct values

Table 2: Ct values obtained for each dilution of the sample control.

Target Name	Sample Name	Reporter	Ct values
<i>Tenacibaculum dicentrarchi</i>	1	FAM	14,12
	2	FAM	17,33
	3	FAM	20,74
	4	FAM	23,8
	5	FAM	26,97
	6	FAM	30,53
	7	FAM	33,82
	8	FAM	-
	Negative Control	FAM	Undetermined

Target Name	Sample Name	Reporter	Ct values
<i>Tenacibaculum maritimum</i>	1	CY5	14,62
	2	CY5	18,29
	3	CY5	21,37
	4	CY5	25,01
	5	CY5	28,57
	6	CY5	32,15
	7	CY5	35,74
	8	CY5	-
	Negative Control	CY5	Undetermined

Certificate validated by:

BQ. Benjamín Catalán - BQ. Rodrigo Iturra Barra

Quality Control Managers