

# Human papillomavirus 42 - Real Time DNA

50 /100 /150 tests (Ready to use kit)

HPV-42 DNA Real Time PCR Kit is a screening assay for a rapid and accurate detection of Human papillomavirus 42.

## Principles of the test:

One Step Bioingentech® PCR Kits provide components for “onestep” real time PCR detection in a convenient format that is compatible with both rapid and standard qPCR cycling conditions.

The One Universal qPCR DNA Master Mix include Bioingentech® all reagents for an optimized qPCR. Specific primers and probe for Human papillomavirus 42 are provided in the kit and can be detected through your Real Time thermal cycler by the 5' nuclease PCR detection method. During PCR amplification forward and reverse primers hybridize to the Human papillomavirus 42 target genomic DNA. A fluorogenic probe is included in the same reaction mixture which consists of a DNA probe labeled with a 5-reporter FAM and 3-quencher which can be detected through green channel.

To confirm extraction of a valid biological template, internal control primers and probe mix are included, which consist of a DNA probe labeled with a 5-reporter HEX and a 3-quencher which hybridize inside a specific housekeeping endogenous target gene. During PCR amplification, the probe is cleaved and the reporter dye and quencher are separated. As a result, a fluorescence increase can be detected on a range of Real Time thermal cyclers through yellow channel.

Our kits also include Positive and Negative Controls which are detailed in the FAQ section.

## Principle and use:

This amplification kit has been manufactured by Bioingentech Ltd. to detect Human papillomavirus 42 via Real Time PCR. This allows absolute quantification or a qualitative assay.

Real Time PCR is based on fluorogenic dyes. Ct values between 0 - 36 should be taken as positive. Values between 36 - 40 Ct should be taken as marginal positive. Ct values above 40 must be considered as negative **(for more details see Table 5)**. This kit needs DNA as a which can be isolated from blood, serum, faeces, respiratory fluid, cerebrospinal fluid, digestive system, tissues, egg yolk, milk, swabs, bacterial cultures, cell lines, among others.

**Table 1. Kit Components:**

| Reactions Tubes                                   | 50 test  | 100 test | 150 test |
|---|----------|----------|----------|
| Universal qPCR Master Mix                         | (1 vial) | (1 vial) | (1 vial) |
| Primer, Probes and Internal Control Universal Mix | (1 vial) | (1 vial) | (1 vial) |
| HPV-42 Positive Control                           | (1 vial) | (1 vial) | (1 vial) |
| HPV-42 Negative Control                           | (1 vial) | (1 vial) | (1 vial) |
| PCR grade Water                                   | (1 vial) | (1 vial) | (1 vial) |

**Table 2. Instrument Compatibility:**

|                         |                             |
|-------------------------|-----------------------------|
| * ABI 7300              | * Mastercycler® ep realplex |
| * ABI 7500FAST          | * LightCycler 2.0 and 480   |
| * ABI 7900              | * Mx3000P QPCR System       |
| * AB Step One           | * Mx3005P QPCR System       |
| * AB Step One Plus      | * RotorGene 3000            |
| * Agilent Mx3005P       | * RotorGene 6000            |
| * CFX96 & CFX384        | * RotorGeneQ                |
| * ExiCycler™ 96         | * SLAN® Real-Time PCR       |
| * iQ5 & MyiQ Cycler     | * Smartcycles II            |
| * Illumina Eco          | * CFX Opus 96               |
| * LightCycler Nano      | * CFX Opus 384              |
| * CFX Opus              | * Quantstudio 3             |
| * Bioer Quant gene 9600 | * Quantstudio 5             |

## Procedure:

Please read through the entire procedure before starting.

## Before Starting

- Pulse-spin each tube in a centrifuge before opening.
- Homogenize the solutions for 5 seconds prior to pipetting
- Use different tips in order to avoid cross contamination.
- Use only sterile, RNase, DNAase, and pyrogen-free tips.

## Step 1

Prepare a Master mix according to the following reaction table 3.

**Table 3. Reaction components for PCR**

| Reaction Tubes                                    | Sample and Internal Control | Positive Control | Negative Control |
|---|-----------------------------|------------------|------------------|
| Universal qPCR Master Mix                         | 7,5 µL                      | 7,5 µL           | 7,5 µL           |
| Primer, Probes and Internal Control Universal Mix | 0,6 µL                      | 0,6 µL           | 0,6 µL           |
| PCR grade Water                                   | 3,4 µL                      | 3,4 µL           | 3,4 µL           |
| DNA Sample  | 3,5 µL                      |                  |                  |
| HPV-42 Positive Control                           |                             | 3,5 µL           |                  |
| HPV-42 Negative Control                           |                             |                  | 3,5 µL           |
| Total Volume                                      | 15 µL                       | 15 µL            | 15 µL            |

## Step 2

Place the tubes in a thermal cycler and perform One Step qPCR according to the program outlined in **Table 4**.

**Table 4. Recommended PCR Cycling table**

| Cycle     | Steps        | Time   | Temp.(°C) |
|-----------|--------------|--------|-----------|
| 1 Cycle   | Denaturation | 3 min  | 95°C      |
| 40 Cycles | Denaturation | 10 seg | 95°C      |
|           | Annealing    | 1 min  | 60°C      |
|           | Hold         | ∞      | 4°C       |

Note1: Collect data in 60 °C

Note 2: Hold Step is optional.

## Interpretation of the test

### 1) Qualitative analysis:

Ct (Threshold cycle) value of each sample can be interpreted as follows.

**Table 5. Ct value result**

| Ct value | Result            |
|----------|-------------------|
| 0 - 36   | Positive          |
| 36 - 40  | Marginal Positive |
| > 40     | Negative          |

1.1.-Ct values over 40 are considered Negative. If the Ct value is in the 0 - 36 range, it must be considered as Positive. This is depending of the sample initial concentration use d for each reaction. Note that the sample real concentration could be altered by the sample purity when it is quantified.

1.2.- For more technical information you must request the **QUALITY CONTROL** for each kits. Also you canrequest more information writing to our email [info@kitpcr.com](mailto:info@kitpcr.com)

## 2) Quantitative analysis:

**Table 6. Preparation of standard curve dilution series. HPV-42 positive control:**

| Average Positive Control Concentration |                     |
|--|---------------------|
| HPV-42                                 | See Quality Control |

| Standar curve | Preparation series a fresh dilution                              | Concentration       | Copy Number         |
|---------------|--|---------------------|---------------------|
| Tube N°1:     | 2uL HPV-42 Positive Control (2 ng/µL) + 18 µL de PCR grade Water | See quality control | See quality control |
| Tube N°2:     | 2uL Tube N°1 + 18 µL de PCR grade Water                          | See quality control | See quality control |
| Tube N°3:     | 2uL Tube N°2 + 18 µL de PCR grade Water                          | See quality control | See quality control |
| Tube N°4:     | 2uL Tube N°3 + 18 µL de PCR grade Water                          | See quality control | See quality control |
| Tube N°5:     | 2uL Tube N°4 + 18 µL de PCR grade Water                          | See quality control | See quality control |
| Tube N°6:     | 2uL Tube N°5 + 18 µL de PCR grade Water                          | See quality control | See quality control |
| Tube N°7:     | 2uL Tube N°6 + 18 µL de PCR grade Water                          | See quality control | See quality control |

2.1. Important Note: Don't forget Homogenize the tubes.

2.2. We will send a Quality Control report for each purchase.

2.3. If you want to obtain less DNA copies you must include a new dilution tube (**Tube N°8**).

2.4. Final DNA copy number will depend of the DNA concentration (**you can see it in Quality Control Report**).

**Table 7. Standard curve set up**

For each of the tubes (tube 1 to tube 7), take 3,5 µl of the final 15 µl mixture.

|   | Tube A | Tube B | Tube C | Tube D | Tube E | Tube F | Tube G |
|---|--------|--------|--------|--------|--------|--------|--------|
| Universal qPCR Master Mix               | 7,5 µL | 7,5 µL | 7,5 µL | 7,5 µL | 7,5 µL | 7,5 µL | 7,5 µL |
| Primer, Probes and Internal Control Mix | 0,6 µL | 0,6 µL | 0,6 µL | 0,6 µL | 0,6 µL | 0,6 µL | 0,6 µL |
| PCR grade Water                         | 3,4 µL | 3,4 µL | 3,4 µL | 3,4 µL | 3,4 µL | 3,4 µL | 3,4 µL |
| Tube N° 1 (Positive Control)            | 3,5 µL |        |        |        |        |        |        |
| Tube N° 2                               |        | 3,5 µL |        |        |        |        |        |
| Tube N° 3                               |        |        | 3,5 µL |        |        |        |        |
| Tube N° 4                               |        |        |        | 3,5 µL |        |        |        |
| Tube N° 5                               |        |        |        |        | 3,5 µL |        |        |
| Tube N° 6                               |        |        |        |        |        | 3,5 µL |        |
| Tube N° 7                               |        |        |        |        |        |        | 3,5 µL |
| Total Volume                            | 15 µL  | 15 µL  | 15 µL  | 15 µL  | 15 µL  | 15 µL  | 15 µL  |

Take table N°6

2.5.- Assess the Ct value when the amplification curve of Standard tubes 1, 2, 3, 4, 5, 6, 7 passes the threshold line. Four tubes are sufficient for a standard curve (tubes 1 to 4).  
 2.6.- Calculate the quantitative value in order to compare the Ct value of the unknown samples with the standard curve.  
 2.7.- When visualizing the result for the Positive Control in the Real Time thermal cycler, there will be only one amplification curve for the target gene. It will not have an amplification curve for the Internal Control.

### 3) Test validation:

3.1.- Each Ct value standard should be as follows. Standard 1 < Standard 2 < Standard 3 < Standard 4 < Standard 5 < Standard 6 < Standard 7.  
 3.2.- R-value of standard curve should be 0.900 - 0.999. R-value represent how well the experimental data fit the regression line. A significant difference in observed Ct values between replicates will lower the R-value.  
 3.3.- The standard curve slope result should be all negative.  
 3.4.- The desired amplification efficiencies vary from 90% to 110%. The theoretical maximum of 100% indicates that the polymerase enzyme is functioning at its maximum capacity. Low reaction efficiencies may be caused by poor primer design or by suboptimal reaction conditions. Reaction efficiencies >110 may indicate pipetting error in your serial dilutions or coamplification of non specific products, such as primer-dimers.

### Visual explanation FAQ:

#### 1.-Positive control:

The Positive control assay uses a FAM dye and should be detected through the green channel of your Real Time thermal cycler (see Table 8).

For copy number determination and as a positive control for the PCR setup, the kit contains a plasmidial positive control template. This can be used to generate a standard curve of Human papillomavirus 42 copy number / Ct value. Alternatively, the positive control can be used at a single dilution if full quantitative analysis of the sample is not required. Each time the kit is used, at least one positive control reaction must be included in the run. We recommend running a positive control for each 12 samples.

A positive result indicates that the primer and probes for the target gene worked properly in that particular experimental scenario. If a negative result is obtained the test results should be invalid and must be repeated (see Table 11). Seal all other samples and negative controls before pipetting the positive control into its well.

### 2.-Internal Control:

The internal control is included in the Primer, Probes and Internal Control Mix alongside the target pathogen detection. The internal control assay uses a HEX dye and should be detected through the Yellow channel of your Real Time PCR instrument and should result in a Ct value of 28 (+/-5), depending on the level of sample dilution and concentration. The refore, a positive result through the yellow channel indicates that PCR conditions are suitable for detection of the target pathogen gene.

If a negative result is obtained through the green channel, the results should be analyzed following the **Table 11 data**.

### 3.-Negative control:


To confirm absence of contamination a negative control reaction should be included every time the kit is used. **Particularly, due to amount of this reagent, you should run a negative control for each 12 samples.** In this instance the Human papillomavirus 42 Negative Control should be used in place of template. A negative result indicates that the reagents have not become contaminated. If a positive result and Ct value less than 36 is obtained, the results should be analyzed and check if a correct amplification curve was obtained. When you obtain a clear amplification curve you should consider repeat your assay due to probably the sample was contaminated (see Table 11).

\* Remember: We recommend running a Positive control and Negative control for each 12 samples. For reaction mix you must use Universal qPCR Master Mix.

### Table 8. Fluorogenic probes, Channels and Dyes

\* In order to setup your software to read the probes, (see Table 8 and 9) and run the corresponding channel:

| Channel | Source   | Detector | Dyes  |
|---------|----------|----------|---|
| Green   | 470 (nm) | 520 (nm) | FAM, Sybr green1, Fluorescein, Eva green, Alexa fluor 488 |
| Yellow  | 530 (nm) | 550 (nm) | Joe, Vic, HEX, Tet, Cal Fluor gold 540, YaKima Yellow     |
| Orange  | 585 (nm) | 610 (nm) | Rox, Cal Fluor Red 610, Cy3.5, Texas Red, Alexa Fluor 568 |
| Red     | 625 (nm) | 660 (nm) | Cy5, Quasar 670, Lightcycler, Red 640, Alexa Fluor 633    |
| Crimson | 680 (nm) | 710 (nm) | Quasar 705, Lightcycler Red 705, Alexa Fluor 680          |

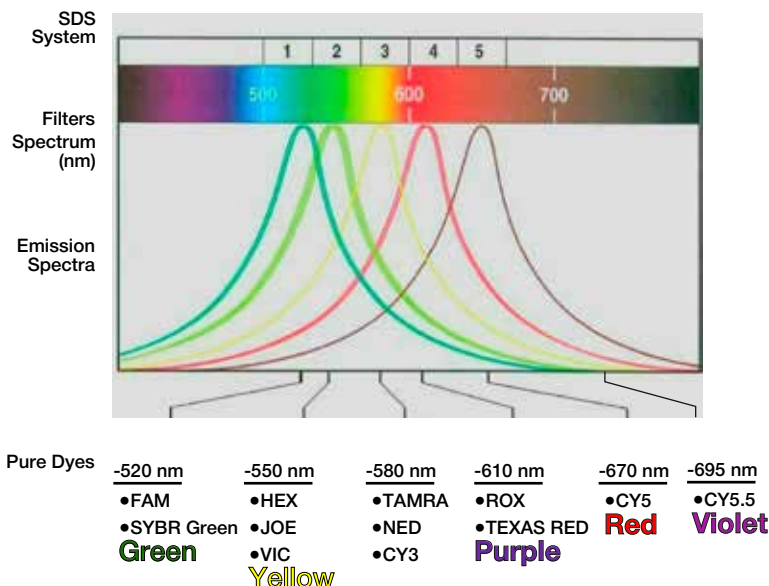

**Table 9. Similarity of our fluorophores with HEX and FAM dyes.**

|                      | Channel | Dyes |
|----------------------|---------|------|
| 1.- Positive Control | Green   | FAM  |
| 2.- Internal Control | Yellow  | HEX  |
| 3.- Sample Target    | Green   | FAM  |

### Important Note:

Request a Certificate of Analysis when you purchase one of our PCR Kits.

Is important mentioned that we can develop special request for other pathogens or multiplex detection according client's specifications.

We strongly recommend don't use or combine our products with reagents from another kits or unknown provenance. We can't assure good result if incompatibility problems occur.

### Temperature

All our reagents are made through protein engineering and are stable at room temperature, the label temperature is just a recommendation after the product is open.

**Table 10. Store Temperature Kits**

|  |                      |                       |
|--|----------------------|-----------------------|
|  | Store temperature    | The label temperature |
|  | Shipping temperature | At room temperature   |

**Table 11. Interpretation of Results**

| Sample | I.C | N.C | P.C | Result                              |
|--------|-----|-----|-----|-------------------------------------|
| +      | +   | -   | +   | POSITIVE                            |
| +      | -   | -   | +   | POSITIVE                            |
| +      | +   | +   | +   | Check Ct and Consider repeat assay* |
| +      | -   | +   | +   | Check Ct and Consider repeat assay* |
| +      | -   | -   | -   | NEGATIVE                            |
| +      | +   | -   | -   | NEGATIVE                            |
| -      | +/- | +/- | +/- | NEGATIVE                            |

\* Sometimes amplification curves for Negative or Internal control with Ct < 30 can be seen, however it might not be a Positive result. If the amplification curve is not sigmoid, you should consider it as a Negative result.

**Table 12. Other Products**

| Products   |
|--|
| OneHumqPCR-realtime™ HPV-42<br>Human papillomavirus 42 Real Time<br>50 tests (Ready to use kit) / Cat. No: Oneq-H471-50D   |
| OneHumqPCR-realtime™ HPV-42<br>Human papillomavirus 42 Real Time<br>100 tests (Ready to use kit) / Cat. No: Oneq-H471-100D |
| OneHumqPCR-realtime™ HPV-42<br>Human papillomavirus 42 Real Time<br>150 tests (Ready to use kit) / Cat. No: Oneq-H471-150D |

**Table 13. Products Specifications**

|                          |  |
|--------------------------|--|
| Technology               | 5' nuclease probe based real time PCR assay  |
| Type of nucleic acid Kit | DNA  |
| Kit storage              | Shipped at room temperature, the label temperature is just a recommendation after the product is open.   |
| Detection Limit          | <b>See Quality Control file. Request it!</b>   |
| Sensitivity & Specify    | Ct value between 0 – 36 should be taken positive. Value between 36-40 Ct should be taken as marginal positive. Ct values above 40 must be considered as negative.      |
| Controls included        | Internal control, Positive control and Negative control included.  |
| Channels                 | <b>FAM Green channel detect pathogen aplicons</b><br><b>HEX Yellow channel detect internal control amplicons.</b><br><b>FAM Green channel detect Positive Control.</b> |